



*Bright Ideas Brought To Light*



## The brightest LED laryngoscope handle

Over 6X brighter than handles using halogen technology.

Fits most fiber optic laryngoscope blades on the market.

Enhanced viewing ability for difficult intubations.



**5 YEAR**  
WARRANTY



## The GreenStar™ laryngoscope handle provides the following benefits:

- **Longer battery run time.** Over 20 hours of operation (4,000 - 5,000 intubations!) from one pair of alkaline batteries.
- **Light intensity doesn't dim** as battery charge diminishes
- Save buying hundreds of lamps during the lifetime of the device
- High illumination levels **extend the usable life** of fiber-optic blade bundles
- **High light output allows enhanced viewing** ability for difficult intubations
- Handles available in **"C" (Medium) size, "AA (Penlight)** and **Stubby**
- **Conforms to ISO 7376-3 standards** for compatibility with all major name brands of green standard fiberoptic laryngoscope blades
- **Tested & certified** to CSA 22. 2 No. 601,1-M & UL 60606-1 Standards
- **304 stainless steel** construction
- **Accessories** see Stellar Series™ fiber optic laryngoscope blades

5 YEAR  
WARRANTY



# Here's how we've put the Green in GreenStar™ to deliver performance and savings:

## Long-life LED

With an average rated life of over **50,000 hours**, there are no more costly halogen lamps to purchase.

## Increased Useable Battery Life – Over 8X Longer\*†

Despite the GreenStar™ laryngoscope handle being the **brightest FO laryngoscope handle available**, it uses less power than a conventional halogen handle. This means you get over 4000 intubations from a single pair of batteries, saving replacement costs and reducing hazardous waste in land-fill sites!


## High Illumination Levels – 6X Brighter\*

This not only benefits the end-users by increasing their viewing ability, the high light output allows blades or bundles that would once have been discarded due to degraded throughput to be successfully used for many more procedures.

\*Over conventional halogen laryngoscope handles

† When comparing the useable intensity of light needed for functional performance of purpose

## GreenStar™ laryngoscope handle performance vs other leading handle

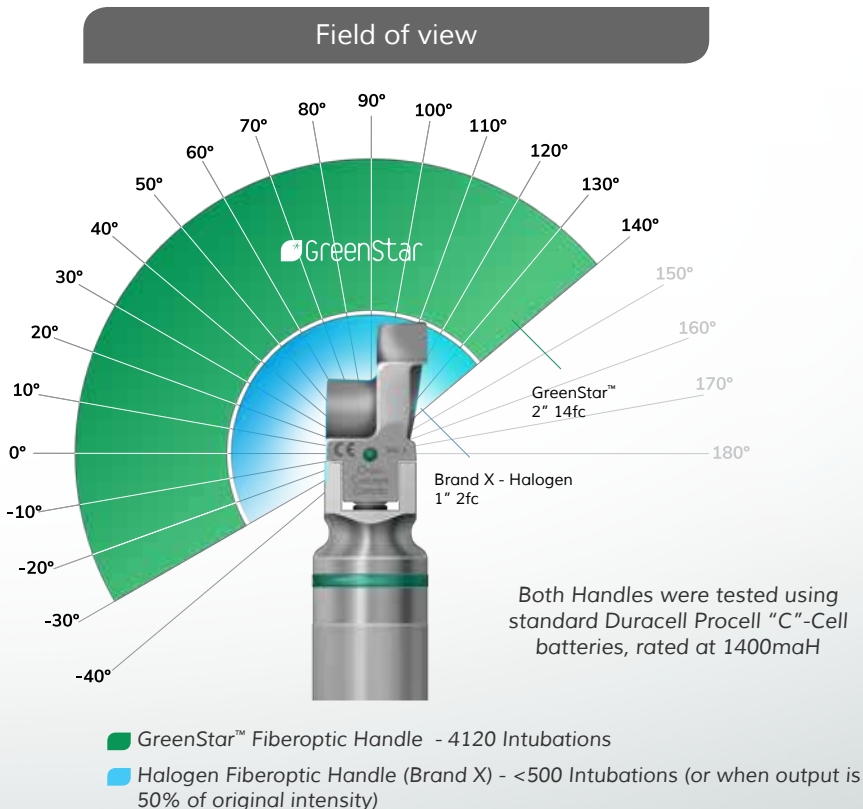

## Upgrade to a brighter, lighter, superior LED Laryngoscope Handle for astounding performance.

The GreenStar laryngoscope family offers the ultimate value package in multiple ways.

- **The brightest LED fiber-optic laryngoscope handle on the market – over 6X the illumination levels** of conventional halogen systems
- Steady light output with **no dimming as batteries age**
- True-colour LED's for accurate visualization
- **Over 20 hours of battery life** using only standard alkaline batteries
- Conforms to ISO 7376-3 Standard (green mark)
- Precisely engineered to offer **perfect compatibility** to other manufacturers products, including Welch Allyn, Heine, Storz and many others
- Surgical grade, 304 Stainless steel construction for **lighter, stronger and more durable** functionality

Distributed by:

 **ORION**<sup>™</sup> [orionconceptsinc.com](http://orionconceptsinc.com)  
C O N C E P T S

Manufacturer of Stellar Series<sup>™</sup> fiber optic laryngoscope blades.